




Infinity Mini Brackets

Features & Benefits

- The Infinity Mini Bracket is a low profile, two-piece bracket with a laser welded mesh ensuring exceptional bonding strength and less breakages and ultimate patient comfort.
- Brackets feature offset bases on the upper and lower pre-molars for optimal bracket positioning.
- Available in .018" Roth & .022" M.B.T / .018" M.B.T & .022" Roth available on request

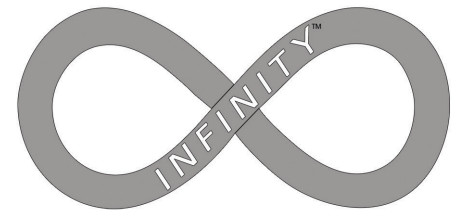


infinity brackets | infinite possibilities

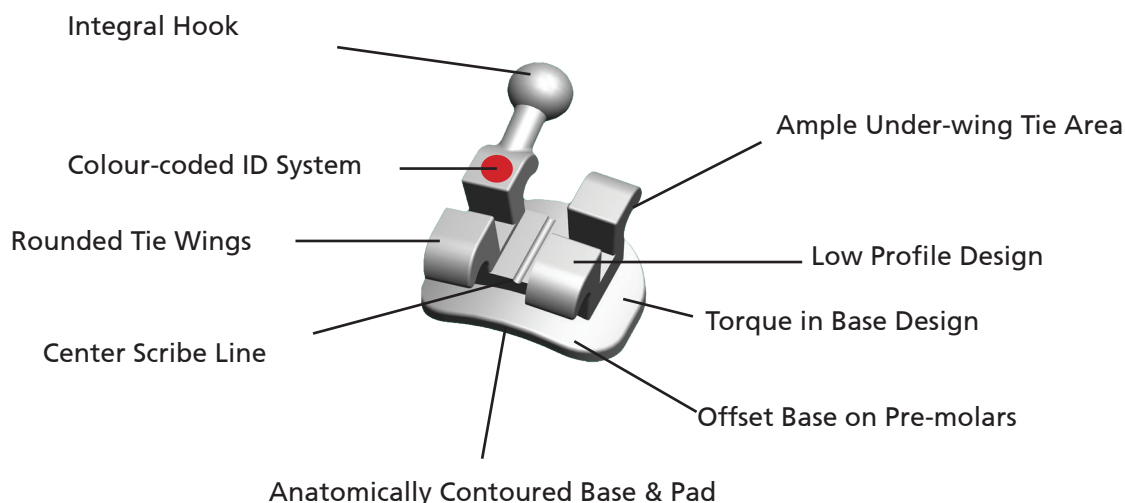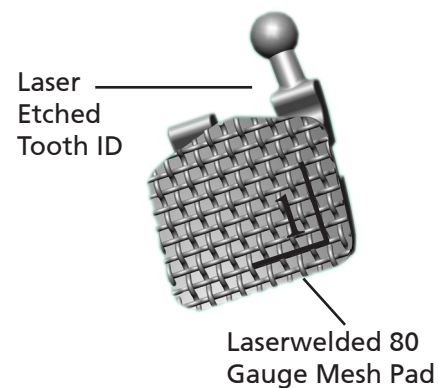

Mini MBT Brackets

in/out			0.7	0.7	0.7	1.08	0.7	0.7	1.08	0.7	0.7	0.7
torque			-7°	-7°	0°	10°	17°	17°	10°	0°	-7°	-7°
tip			0°	0°	8°	8°	4°	4°	8°	8°	0°	0°
width			2.8	2.8	2.8	2.8	3.1	3.1	2.8	2.8	2.8	2.8
Tooth	18/22	w/ hook	U5RHK	U4RHK	U3RHK	U2R	U1R	U1L	U2L	U3LHK	U4LHK	U5LHK
		w/o hook	U5R	U4R	U3R					U3L	U4L	U5L

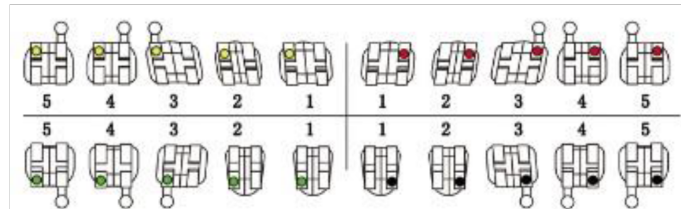


in/out			0.7	0.7	0.7	0.9	0.9	0.9	0.9	0.7	0.7	0.7
torque			-17°	-12°	-0°	-6°	-6°	-6°	-6°	-0°	-12°	-17°
tip			2°	2°	3°	0°	0°	0°	0°	3°	2°	2°
width			2.8	2.8	2.8	2.2	2.2	2.2	2.2	2.8	2.8	2.8
Tooth	18/22	w/ hook	L5RHK	L4RHK	L3RHK	L2R	L1R	L1L	L2L	L3LHK	L4LHK	L5LHK
		w/o hook	L5R	L4R	L3R					L3L	L4L	L5L

Size	Package	Hooks	
18/22	With laser mark set	3 with hook	345 with hook

Please note: Minimum quantity order for single brackets is 5 per tooth

Mini Roth Brackets

in/out			0.68	0.68	0.66	1.06	0.68	0.68	1.06	0.66	0.68	0.68
torque			-7°	-7°	-2°	8°	12°	12°	8°	-2°	-7°	-7°
tip			0°	0°	11°	9°	5°	5°	9°	11°	0°	0°
width			2.8	2.8	3.2	2.8	3.2	3.2	2.8	3.2	2.8	2.8
Tooth	18/22	w/hook	U5RHK	U4RHK	U3RHK	U2R	U1R	U1L	U2L	U3LHK	U4LHK	U5LHK
		w/o hook	U5R	U4R	U3R					U3L	U4L	U5L

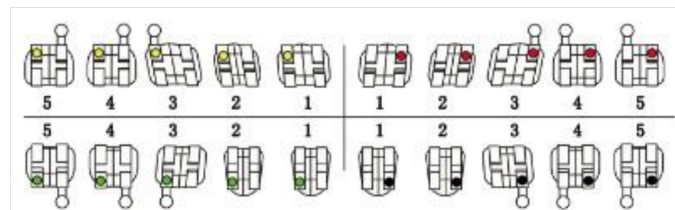


in/out			0.68	0.72	0.66	1.06	1.06	1.06	1.06	0.68	0.72	0.68
torque			-22°	-17°	-11°	-1°	-1°	-1°	-1°	-11°	-17°	-22°
tip			0°	0°	5°	0°	0°	0°	0°	5°	0°	0°
width			2.8	2.8	2.8	2.5	2.5	2.5	2.5	2.8	2.8	2.8
Tooth	18/22	w/hook	L5RHK	L4RHK	L3RHK	L2R	L1R	L1L	L2L	L3LHK	L4LHK	L5LHK
		w/o hook	L5R	L4R	L3R					L3L	L4L	L5L


Size	Package	Hooks	
18/22	With laser mark set	3 with Hook	345 with hook

Please note: Minimum quantity order for single brackets is 5 per tooth

Infinity Plus Brackets

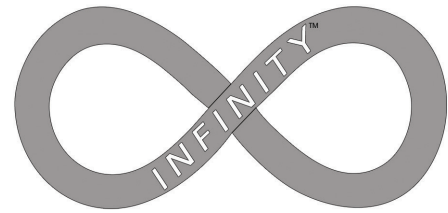

infinity brackets | infinite possibilities

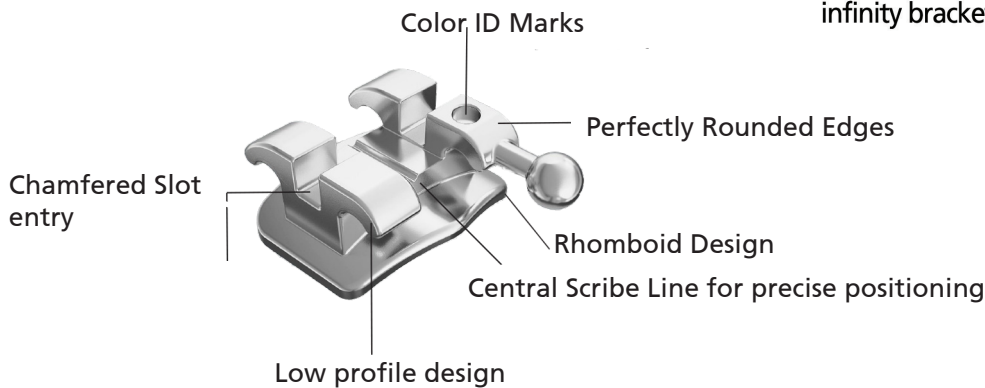

Sandblasted base provides excellent adaptation and adhesion



Features & Benefits

- Infinity Plus low-profile brackets have all the advantages of our Infinity Mini Brackets but also have a mesh base (rather than a cast MIM base) which has also been sandblasted thereby further increasing retention.
- They have a chamfered slot and perfectly round edges for patient comfort.
- Manufactured from the best quality 17-4 Stainless Steel and MIM (Metal Injection Mould) technology.
- Available in mini or standard size.
- Available in .018" / .022" Roth & M.B.T.