




© Photo: J Kurtz-Hoffmann

# Photo-Mirror-Demister


## The Mist-Free Mirror For Dental Photography



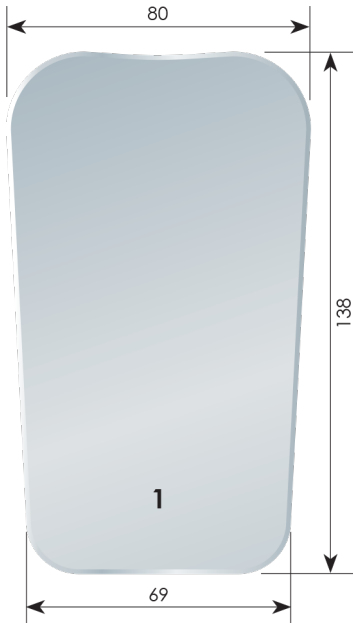


**FSFDM01 R1495 COD**

*Stainless Steel Photo Mirror, Front Surface, double side mirror surface  
Occlusal large, Adults*

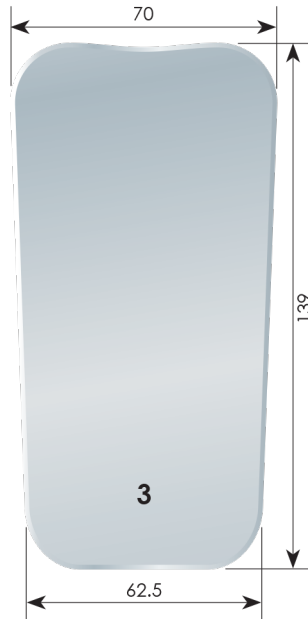


**FSFDM03 R1495 COD**

*Stainless Steel Photo Mirror, Front Surface, double side mirror surface  
Occlusal standard, Adults*

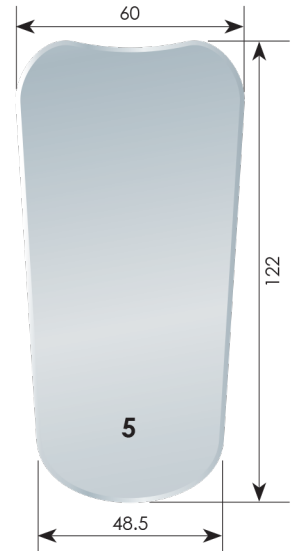


**FSFDM05 R1495 COD**

*Stainless Steel Photo Mirror, Front Surface, double side mirror surface  
Occlusal small for Children*

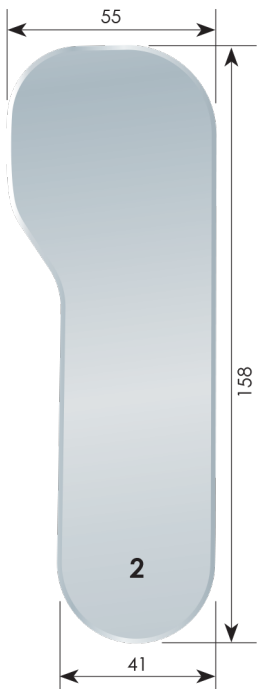


**FSFDM02 R1495 COD**

*Stainless Steel Photo Mirror, Front Surface, double side mirror surface  
Lateral, Oral, Adults*

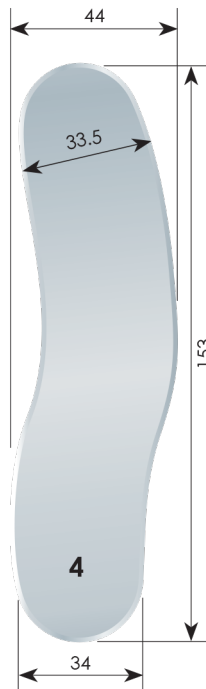

**FSFDM04 R1495 COD**

*Stainless Steel Photo Mirror, Front Surface, double side mirror surface  
Lateral, for Children*

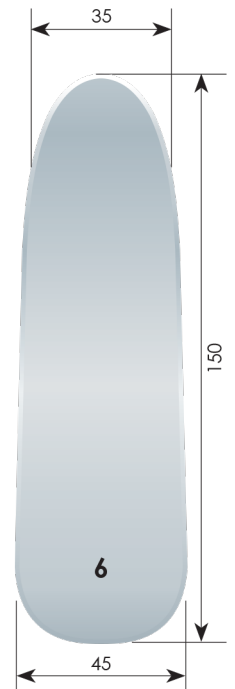

**FSFDM06 R1495 COD**

*Stainless Steel Photo Mirror, Front Surface, double side mirror surface  
Lateral, Universal*

## PMDECO R7500 COD

ECONOMY kit consisting of:

- Mirror holder / main unit
- Lithium battery and charger
- Operating instructions
- 1 Photo Mirror of your choice

### BULK DEAL

Buy 2 or more  
and only pay  
**R6800 COD**



## PMDECO+ R8900 COD

ECONOMY PLUS kit consisting of:

- Mirror holder / main unit
- Lithium battery and charger
- Operating instructions
- 1 Photo Mirror of your choice
- 1 Vestibule retractor for CRLA
- 2 Cheek retractors Columbia CRC

### BULK DEAL

Buy 2 or more  
and only pay  
**R7900 COD**



## PMDPRO R21900 COD

PROFESSIONAL kit consisting of:

- Suitcase for photo equipment
- Mirror holder / main unit
- Lithium battery and charger
- Operating instructions
- Photo Mirror Tray PMDK03
- 4 Photo Mirrors of your choice
- 1 Vestibule retractor for CRLA
- 2 Cheek retractors Columbia CRC

### BULK DEAL

Buy 2 or more  
and only pay  
**R19750 COD**



## PMDEVE R16500 COD

EVERY kit consisting of:

- Mirror holder / main unit
- Lithium battery and charger
- Operating instructions
- Photo Mirror Tray PMDK03
- 4 Photo Mirrors of your choice

### BULK DEAL

Buy 2 or more  
and only pay  
**R14800 COD**



## PMDVAR R14900 COD

VARIO kit consisting of:

- Mirror holder / main unit
- Lithium battery and charger
- Operating instructions
- Photo Mirror Tray PMDK03
- 3 Photo Mirrors of your choice

### BULK DEAL

Buy 2 or more  
and only pay  
**R13500 COD**



## PMDVAR+ R16200 COD

VARIO PLUS kit consisting of:

- Mirror holder / main unit
- Lithium battery and charger
- Operating instructions
- Photo Mirror Tray PMDK03
- 3 Photo Mirrors of your choice
- 1 Vestibule retractor for CRLA
- 2 Cheek retractors Columbia CRC

### BULK DEAL

Buy 2 or more  
and only pay  
**R14500 COD**



# MIST-FREE PHOTO MIRRORS

## *Features and Benefits*

### **Mist free surface of photo mirror due to the air ventilation**

- Air volume can be increased or decreased according to patient requirement and/or doctor preference via +/- keypad.
- The air volume level will be indicated on the keypad via LED's.
- The last position of air volume level (after 30 seconds use) will be memorized for the next use.

### **Angulation of the Photo Demister Handle**

- More comfortable handling.
- Better access for the digital camera.
- Can be easily held by patient (doctor can take photos without assistant).

### **Modified mirror clamp fixation system, made of stainless steel**

- Safe holding of the mirrors.
- Easy One-Screw-Fixation of the mirror.
- Photo Mirrors can be fixed in the PMD at every side,

### **New material of PMD-Corpus**

- Contains silver ions - antibacterial influence.
- Can be completely disinfected with a disinfection cloth.
- No discolouration/ changing of colour even after extended use.
- Light weight.
- More robust construction.

### **The LED light**

- Helps focusing of digital camera.
- No blinding of patient eyes.
- Dispersive light beam.

### **Changing of battery**

- Uncomplicated and very easy via standard Micro-USB connector.
- Red LED indication of chargin progress.
- Battery overload protection.
- Now available with EU & USA sockets.
- Easy battery replacement - standard battery pack, used in many digital cameras of known brands.
- Easy charger replacement - charger, used in many smartphones, tablets etc. of known brands.

### **Stainless steel mirrors**

- Better thermal conductivity than by glass mirrors - very quick adoption of patient body temperature.
- More robust and longer life than glass mirrors.
- Mirror surface on both sides.
- Easy to clean and sterilise (also ultrasonic cleaning).

### **Diverse accessories for more comfort**

- Sterilisation cassette for photo mirrors - comfortable sterilisation and storage.
- Retractors made of stainless steel - very long life.

