

MIST-FREE PHOTO MIRRORS

Features and Benefits

Mist free surface of photo mirror due to air stream

- Air volume can be increased or decreased according to patient requirement and/or doctor preference via +/- keypad.
- The air volume level will be indicated on the keypad via LED's.
- The last position of air volume level (after 30 seconds use) will be memorized for the next use.

Angulation of the Photo Demister Handle

- Allows more ergonomic handling.
- Better access for the digital camera.
- Can be easily held by patient (Doctor can take photos without an assistant).

Modified mirror clamp fixation system, made out of stainless steel

- Secure holding of the mirrors.
- Easy One-Screw-Fixation of the mirror.
- Photo Mirrors can be fixed in the PMD on both sides.

New material of Photo Mirror Demister Cover

- Contains silver ions - antibacterial influence.
- Can be completely disinfected with a disinfectant wipe.
- Will not discolour, even after extended disinfection and use.
- Light weight.
- More robust construction.

The LED light

- Helps with focusing of the digital camera.
- Will not blind patients eyes.
- Light beam completely illuminates oral cavity.

Charging of battery

- Very easy via standard Micro-USB connector.
- Red LED indication of charging progress.
- Battery overload protection.
- Easy battery and charger replacement - available to purchase from Orthoshop.

Stainless steel mirrors

- Better thermal conductivity than glass mirrors – therefore less misting up.
- More robust and longer life than glass mirrors.
- Mirror surface is on both sides.
- Easy to clean and sterilise (also ultrasonic cleaning).

Diverse accessories for more comfort

- Sterilisation cassette for photo mirrors - comfortable sterilisation and storage.
- Retractors made of stainless steel - very long life.



EXCLUSIVE TO ORTHOSHOP



© Photo: J Kurtz-Hoffmann

Photo-Mirror-Demister The Mist-Free Mirror For Dental Photography



+27 (0)21 976 5403



sales@orthoshop.co.za



+27 (0)82 907 9697



www.orthoshop.co.za



orthoshop sales






OrthoshopInternational



Jakobi videos on our YouTube channel: search Orthoshop International

PHOTO MIRRORS

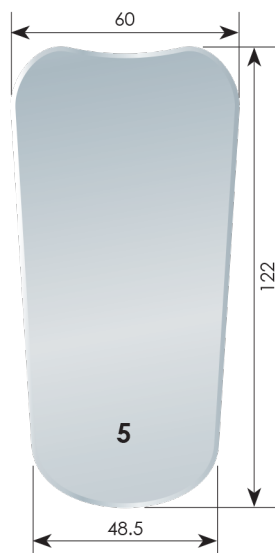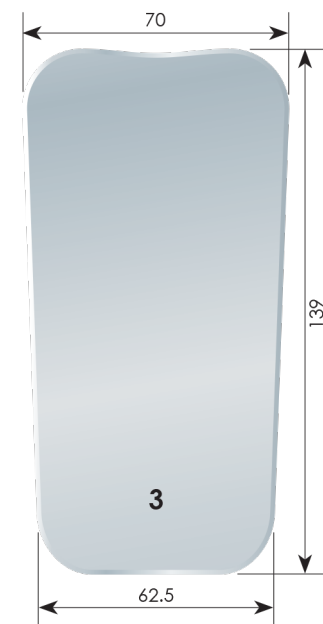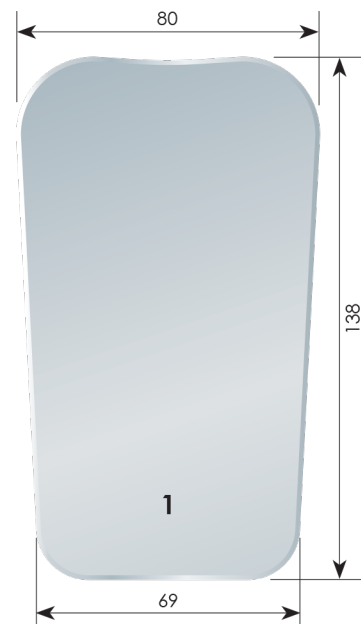

FSFDM01




Stainless Steel Photo Mirror, Front Surface, double side mirror surface
Occlusal large, Adults

FSFDM03

Stainless Steel Photo Mirror, Front Surface, double side mirror surface
Occlusal standard, Adults

FSFDM05

Stainless Steel Photo Mirror, Front Surface, double side mirror surface
Occlusal small for Children

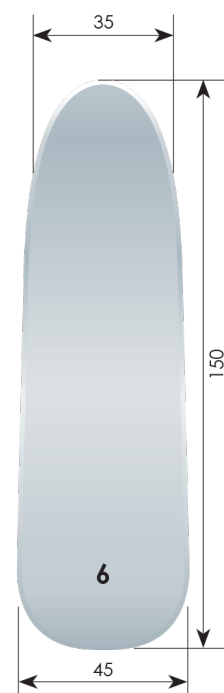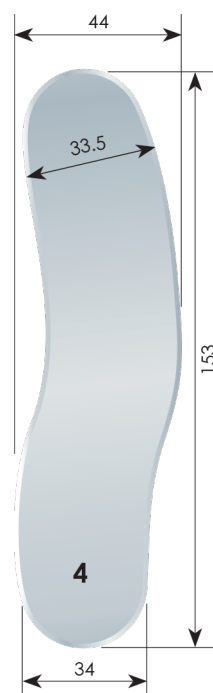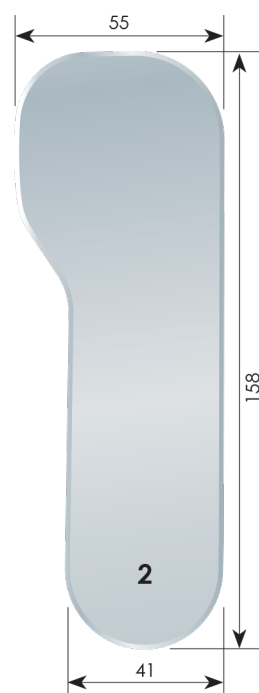

FSFDM02

Stainless Steel Photo Mirror, Front Surface, double side mirror surface
Lateral, Oral, Adults

FSFDM04

Stainless Steel Photo Mirror, Front Surface, double side mirror surface
Lateral, for Children

FSFDM06

Stainless Steel Photo Mirror, Front Surface, double side mirror surface
Lateral, Universal

* Measurements in mm

PHOTO MIRROR KITS

PMDECO

ECONOMY kit consisting of:

- Mirror holder / main unit
- Lithium battery and charger
- Operating instructions
- 1 Photo Mirror of your choice



PMDECO+

ECONOMY PLUS kit consisting of:

- Mirror holder / main unit
- Lithium battery and charger
- Operating instructions
- 1 Photo Mirror of your choice
- 1 Vestibule retractor for CRLA
- 2 Cheek retractors Columbia CRC



PMDPRO

PROFESSIONAL kit consisting of:


- Suitcase for photo equipment
- Mirror holder / main unit
- Lithium battery and charger
- Operating instructions
- Photo Mirror Tray PMDK03
- 4 Photo Mirrors of your choice
- 1 Vestibule retractor for CRLA
- 2 Cheek retractors Columbia CRC



PMDEVE

EVERY kit consisting of:

- Mirror holder / main unit
- Lithium battery and charger
- Operating instructions
- Photo Mirror Tray PMDK03
- 4 Photo Mirrors of your choice

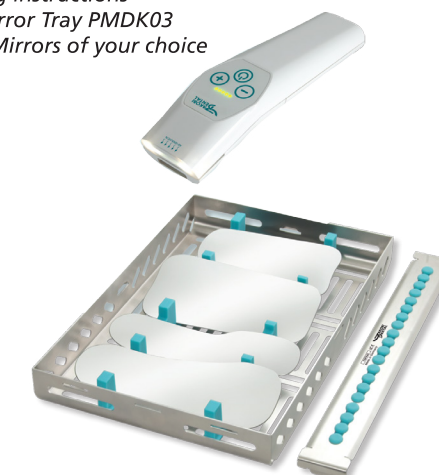

PMDVAR

VARIO kit consisting of:

- Mirror holder / main unit
- Lithium battery and charger
- Operating instructions
- Photo Mirror Tray PMDK03
- 3 Photo Mirrors of your choice



PMDVAR+

VARIO PLUS kit consisting of:

- Mirror holder / main unit
- Lithium battery and charger
- Operating instructions
- Photo Mirror Tray PMDK03
- 3 Photo Mirrors of your choice
- 1 Vestibule retractor for CRLA
- 2 Cheek retractors Columbia CRC

